
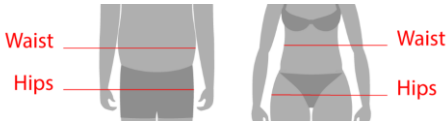




Skin and Wound Product Information Sheet

Knit Pants	
Classification	Continence Containment: Knit Pant
Key Points	<ul style="list-style-type: none"> To used as part of a 2-piece continence containment system; holds the urinary continence pad in place. Pants are made with solid knit 2-way stretch polyester with a spandex waist and has a seamless design which eliminate seams pressing into the skin. Latex-free. Pants can be machine washed multiple times. Color-coded waistbands make size identification easier.
Indications	<ul style="list-style-type: none"> Use in conjunction with a pad insert for clients who have urinary continence and/or to support toilet retraining.
Precautions	<ul style="list-style-type: none"> Correct size is required to hold the pad in place and prevent skin damage in groins and waist area.
Contraindications	<ul style="list-style-type: none"> Do not double pad. Pads are not designed to allow excess urine to flow into the pad beneath. The extra layer can contribute to heat/humidity within the pad and pressure leading to Incontinence Associated Dermatitis (IAD) and/or Pressure Injury.
Formats & Sizes	<p>Sizes</p> <ul style="list-style-type: none"> ○ Medium/Large: blue/grey stripe Waist size 51-152cm (20-61in) ○ XL: green stripe Waist size 114-178cm (45-70in) ○ 2XL: purple stripe Waist size 127-191cm (51-76in) ○ 3XL: red stripe Waist size 162-213cm (65-85in) <p>There is no Small size.</p> <div style="text-align: right; margin-top: 10px;">  </div>
Application Directions	Rationale
<p>Measure client's waist to ensure correct size of pant.</p> <div style="text-align: center; margin: 10px 0;">  </div> <p>Position urinary continence pad in the pant. Ensure that the pad's adhesive strip (if present) sticks to pants.</p> <p>Position the pad for males versus females; mid-to-front for males and middle for females. Also consider the client's usual body position; middle if mostly sitting and mid-to-back if in bed.</p>	<p>Must ensure that the pad conforms to client's perineal area to prevent leakage.</p>
Daily Care	
Change pants when soiled. Pants can be hand-washed or machine-washed/tumbled-dry.	
Expected Outcome	
Pants supports the pad comfortably with no leakage of urine.	
For further information, please contact your NSWOC or Wound Clinician.	